

Basic benches with power console

BENCH WITH POWER CONSOLE



6 power sockets + 1 circuit breaker + 1 indicator light.
A PVC protective pad is placed underneath each foot.
Maximum load 1000kg uniformly distributed.

Dims (cm)	Thermal magnetic circuit breaker		30mA residual current circuit breaker	
	Ref.	Ref.	Ref.	Ref.
180x75	FSV-187	FSO-187	FSV-187D	FSO-187D
180x100	FSV-1810	FSO-1810	FSV-1810D	FSO-1810D
200x75	FSV-207	FSO-207	FSV-207D	FSO-207D
200x100	FSV-2010	FSO-2010	FSV-2010D	FSO-2010D

BENCH WITH POWER CONSOLE & EMERGENCY STOP



Power console with: 6 power sockets + 1 circuit breaker
+ 1 indicator light + 1 emergency stop with key (positive safety).
A PVC protective pad is placed underneath each foot.
Maximum load 1000kg uniformly distributed.

Dims (cm)	Thermal magnetic circuit breaker		30mA residual current circuit breaker	
	Ref.	Ref.	Ref.	Ref.
180x75	ESV-187	ESO-187	ESV-187D	ESO-187D
180x100	ESV-1810	ESO-1810	ESV-1810D	ESO-1810D
200x75	ESV-207	ESO-207	ESV-207D	ESO-207D
200x100	ESV-2010	ESO-2010	ESV-2010D	ESO-2010D

INFORMATION

All power sockets are wired (The power cable of the console is not supplied)
Recommendation: use a power cable in 3x2.5mm² (phase + neutral + earth)


VERTICAL CABLE RECEWAY

For the wiring of the power cable
from the aerial wiring system
to the base of the power supply

ref. DEG-1

Model for power cable and
network cable

ref. DEG-2

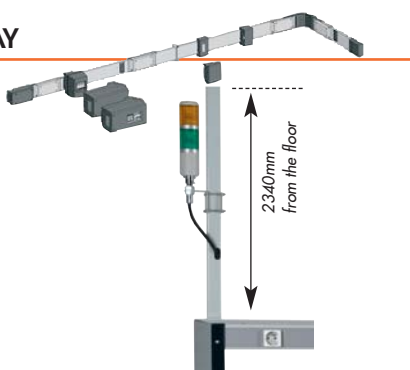

EQUILIBRIUM RANGE

New generation of technical furniture

Quality of steel & finishes



All parts of the metal frame are manufactured from electrogalvanised steel sheet, which is widely used in the car manufacturing industry and makes the sheet extremely resistant to corrosion. The iron/zinc combination forms a natural cell with moisture in the air that prevents the iron from corroding, even if shearing or scratching occurs. This protection is strengthened further by two layers of furnace-baked epoxy paint, which means that it is suitable for use even in tropical or coastal areas.

RAL 2008

Colour Chart

RAL 6018

Products shown in this

RAL 1021

catalogue, especially

RAL 3003

workbenches, can be

RAL 5017

supplied in the colours

RAL 7016

given here. All other RAL

RAL 7035

available upon request.

Quality of stratified top



SHOCK
RESISTANT



ABRASION
RESISTANT



SCRATCH
RESISTANT



CHEMICAL
RESISTANT



EASY
TO CLEAN



Uniformly
distributed

Thickness of 40mm



High quality stratified layer. High heat resistance (up to 300°C without modification in shape, colour). High resistance to scratches

High density particle board
(650 kg/m³) with certified wood
(Environnemental Charter PEFC)

Gluing and
power press

PVC sides - Heat gluing

Stratified layer of stabilization (avoid
bending & torsion, whatever is the climate)

Power sockets available for technical furniture



F TYPE SCHUKO
GERMANY



E TYPE
FRANCE



J TYPE
SWITZERLAND



G TYPE
UNITED KINGDOM



B TYPE
USA